



200



1200

Injection Techniques

Sunekos® 200 is easy to use and is compatible with most injection techniques.

Sunekos® 200 can be injected with:

A

TRADITIONAL
TECHNIQUES
(MICRO-BOLUS OR
RETROGRADE
LINEAR INJECTION)

B

WITH A
PROPRIETARY
TECHNIQUE TO ACHIEVE
THE BEST RESULTS
(IFT: INTERSTITIAL
FLUID TECHNIQUE)

IFT
INTERSTITIAL
FLUID TECHNIQUE

FACE

NECK

DÉCOLLETÉ

HANDS

EYE CIRCLES AND BAGS

FOREHEAD

INNER THIGH



Sunekos[®] 200 Standard Protocol

INDICATION

- Sunekos[®] 200 is the first line product and it is suitable for all types of dermis.
- 4 treatments of Sunekos[®] 200 at 7/10 day intervals each.



MAINTENANCE

Repeat  after 6 months.

DAILY SUPPORT

For a holistic approach add **Nutrakos Drinkable** every day.




ALTERNATIVE

Alternatively MAINTENANCE can be done with 2 treatments every 3 or 4 months, repeating the full protocol with 4 treatments every 12 months.

Sunekos[®] 200 Micro-Bolus injection areas

- Needle 27-30G and 12-13 mm
- Medium-deep Dermis injection
- Micro bolus or retrograde linear injection
- 0.1 ml injection

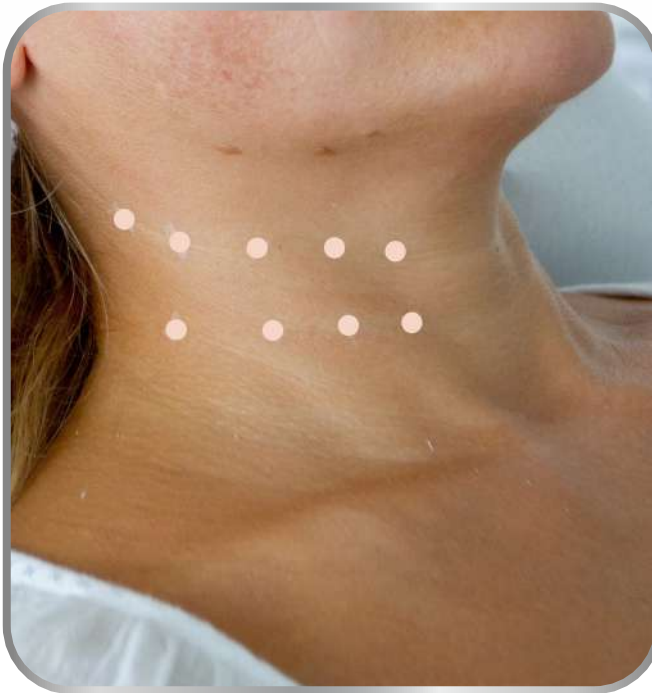

IFT: Interstitial Fluid Technique



- Needles 27-30G and 12-13 mm
- Medium-deep Dermis injection
- 45° Retrograde linear injection
- Series of bolus - Distance of about 1.5 cm
- 0.1 ml of injection
- Injections in the directions of flow
- Massage following the flow

Interstitial Fluid Technique is an injection technique developed by **Professional Dietetics medical team** to get faster results(*) especially in the first phase.



IFT is combining the effect of **2 fundamental actions** for the treatment of the Extra Cellular Matrix in the dermis:

- 1 Medium-deep Dermis injection**
- 2 Lymphatic Drainage**

(*) clinical trial under publication process.

IFT: Interstitial Fluid Technique: Patient Customised Methodology

The **direction** of the sequence of the injection points will allow you to obtain a certain aesthetic effect (**custom infiltrative technique**).

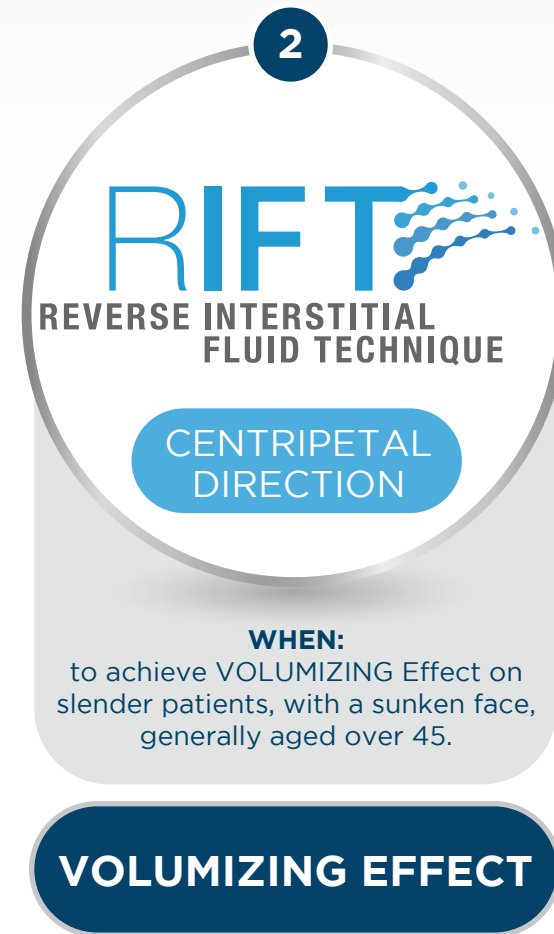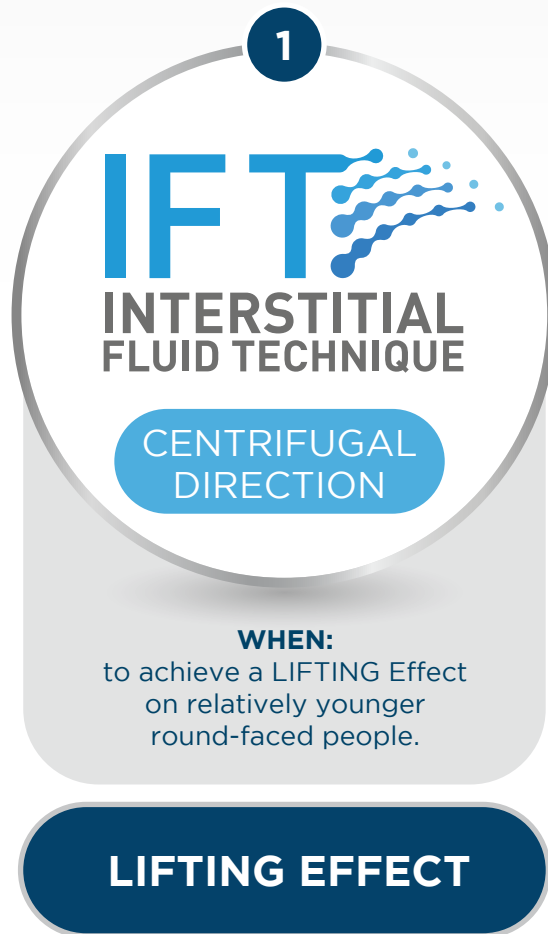

IFT: Interstitial Fluid Technique: Patient Customised Methodology



**INNOVATIVE
PROPRIETARY
TECHNIQUE
LAUNCHED AT AN
INTERNATIONAL
LEVEL**

**APPRECIABLE
VISIBLE RESULTS
CONFIRMED BY
INSTRUMENTAL
EVALUATIONS
WITH *VECTRA***

**METHODS
COMBINED AND
ADAPTABLE TO
THE DIFFERENT
CONFORMATIONS
OF THE FACE**



**CUSTOMIZABLE
TECHNIQUE
BASED ON
DIFFERENT
MECHANISMS OF
SKIN AGING**

IFT: Clinical Trial Results

Interstitial Fluid Technique gives faster results especially in the initial phase:

**2.3 TIMES FASTER
ON CROW'S FEET
SEVERITY RATE
AFTER 7 DAYS**

**2.5 TIMES FASTER
ON WRINKLES
SEVERITY GRADE
AFTER 7 DAYS**

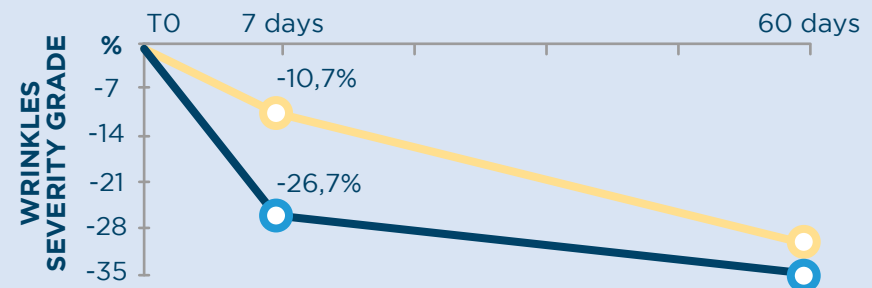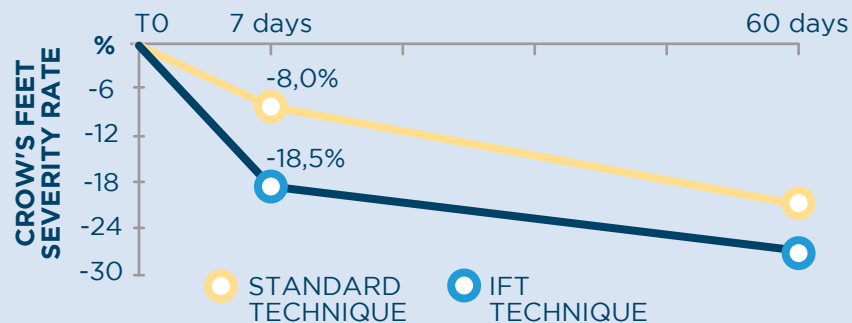


Sunekos[®] 200


IFT
INTERSTITIAL
FLUID TECHNIQUE

Injection areas

- IFT leverages on the flows of interstitial fluid
- Needle 27-30G and 12-13 mm
- Medium-deep Dermis injection
- 45° Retrograde Linear Injection

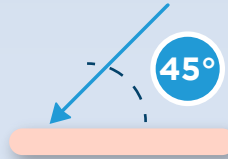

- Injections in the direction of the flow
- 0.1 ml injection
- Distance about 1.5 cm
- Massage following the flow

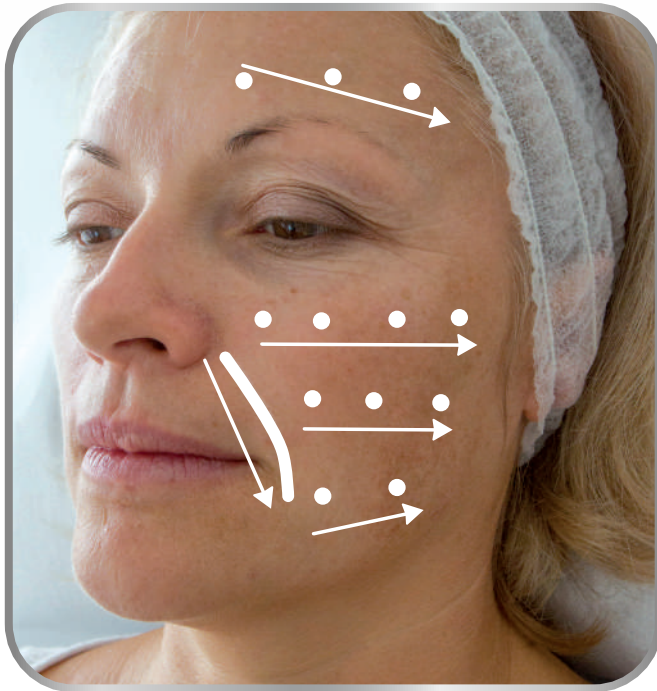

Reverse Interstitial Fluid Technique

**SEVERE WRINKLES
AND SAGGING SKIN**

- **RIFT** is an injection technique with opposite direction to interstitial fluids.
- It is performed in areas of depletion and emptying of volumes.
- The sequence of injections goes from the periphery towards the center.
- Following, a light massage with the fingertip in order to support and stimulate the physiological flow of interstitial fluids, from the periphery towards the center.



RIFT helps to distribute the product on the center-face area and to stimulate the neo-regeneration of Extra Cellular Matrix.

IMMEDIATE VOLUMIZING EFFECT

Severe Wrinkles & Sagging Skin

Sagging Skin and Severe Wrinkles are mainly caused by loss of volume.

To get the best results from **Sunekos®** we have to combine the action of a Medium/High density Hyaluronic Acid with a Low density HA such as **Sunekos® 200**.

The solution is to use **Sunekos® 1200** in combination with **Sunekos® 200** creating a “scaffold” through the “**cushion**” technique.

“cushion” technique



An internal trial comparing 2 groups with sagging skin, the first treated only with S200 and the second with one x S1200 at the beginning as “**cushion**” followed by S200 standard treatment was giving the following results:

Facial Volume
Measurement in CC

+24%

FVLS improvement
(photographic scale)

+33%

Severe Wrinkles & Sagging Skin

Sunekos[®] 200 + Sunekos[®] 1200 Protocol

INDICATION

- This protocol is for **thin/slack faced relatively elder people.**




1 treatment with Sunekos[®] 1200 as a Scaffold followed, on the same day or no more than 1 or 2 days later, by 3 or 4 treatments of **Sunekos[®] 200** at 7/10 day intervals each.



MAINTENANCE

Maintenance with 1 treatment
Sunekos[®] 1200 + 1 treatment **Sunekos[®] 200**

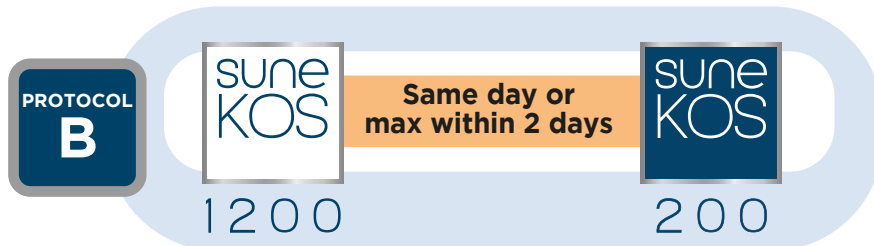

DAILY SUPPORT

For a holistic approach add
Nutrakos Drinkable every day.



Repeat after **4 and 7 months** from first treatment.

Sunekos[®] 200 + Sunekos[®] 1200

Cannula and RIFT Injection areas

SEVERE
WRINKLES

- Cannula 25 or 27G and 50 mm
- Deep Dermis injection
- Fan Technique

STEP
1



sune
KOS

1200


- Needle 27-30G and 12-13 mm
- Medium-deep Dermis injection
- 45° Retrograde Linear Injection
- Needle in direction against the flow
- 0.1 ml injection at 1.5 cm distance
- Massage against the flow

STEP
2



sune
KOS

200

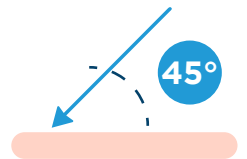

Sunekos® 200 Dark Eye Circles & Bags

THIN SKIN NEAR
THE EYE AREA IS
ONE OF THE
REASONS FOR
**DARK EYE
CIRCLES**



200

IN THIS
CONDITION
YOU SEE IN
TRANSPARENCY
THE TISSUES
BELOW




Sunekos® 200, injected with a specific technique helps to stimulate the neo-regeneration of Extra Cellular Matrix giving thickness to the Skin.

QUICK COVERING EFFECT

Sunekos® 200

Technique Dark Eye Circles & Bags

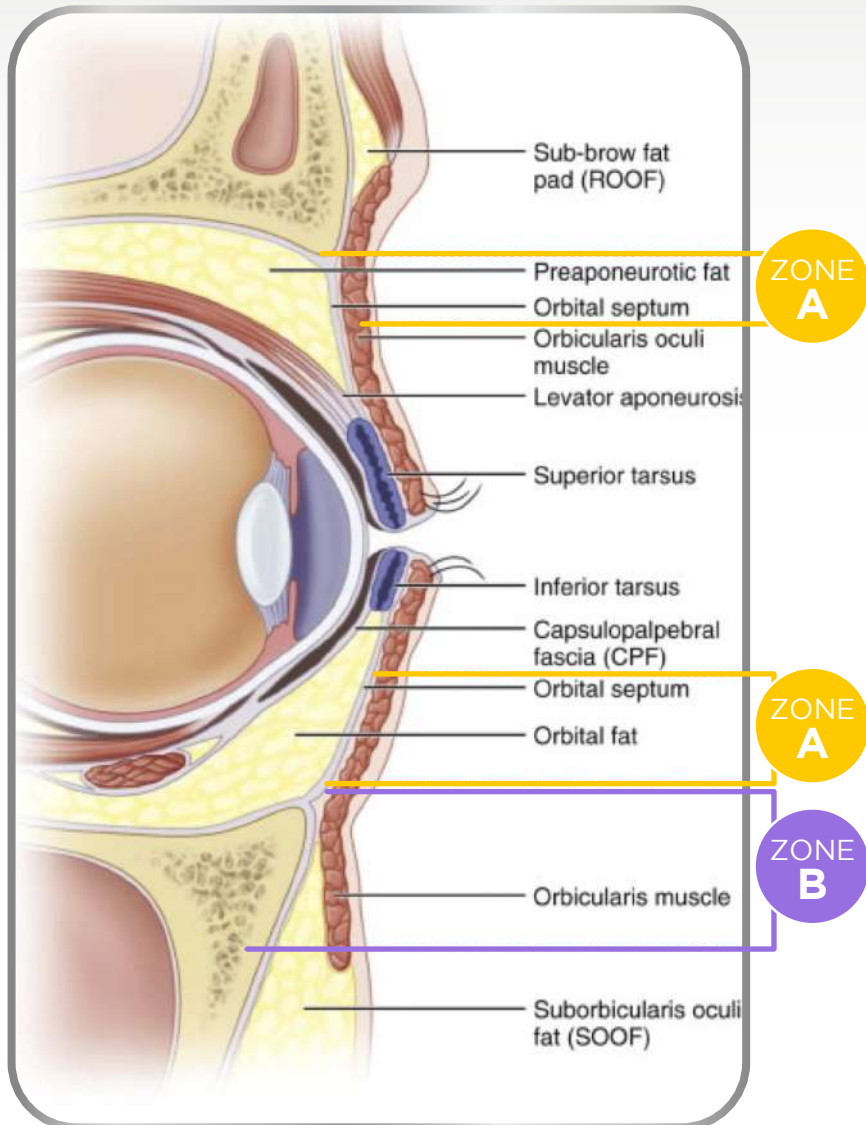

DUE TO THE CAPABILITY TO **PRODUCE NEW ECM**, IS PARTICULARLY EFFECTIVE TO REDUCE DARK EYE CIRCLES.

sune
KOS
200

BEING ABLE TO **REPRODUCE ELASTIN**, IT CAN SIGNIFICANTLY REDUCE EYE BAGS.

- Before injecting in the periocular area, please be mindful of its anatomy.
- Over the orbital septum there is no vascularization.
- Any injection over orbital septum can generate damages.
- NEVER HIT THE ORBITAL SEPTUM. To avoid it please pay attention to:

ZONE A

inject parallel to the skin & very superficially

ZONE B

inject perpendicularly until touching the bone

Clinical trial under publication process.

Sunekos[®] 200

Technique for Dark Eye Circles & Bags

ZONE
A

- Inject parallel to the skin & very superficial.
- Inject on the yellow lines in one 1 or 2 points.

ZONE
B

- Inject perpendicularly till touching the bone.

CAUTION

- Take extra care and avoid hitting veins.
- Carefully check the skin to avoid injecting the small veins, especially **ZONE A**.
- Reduce number of points of injection if the case.

3 SESSIONS

sune
KOS
200

14 days

sune
KOS
200

14 days


sune
KOS
200

+

TO PRESERVE
RESULTS

nutra
KOS

Gel Eye Contour

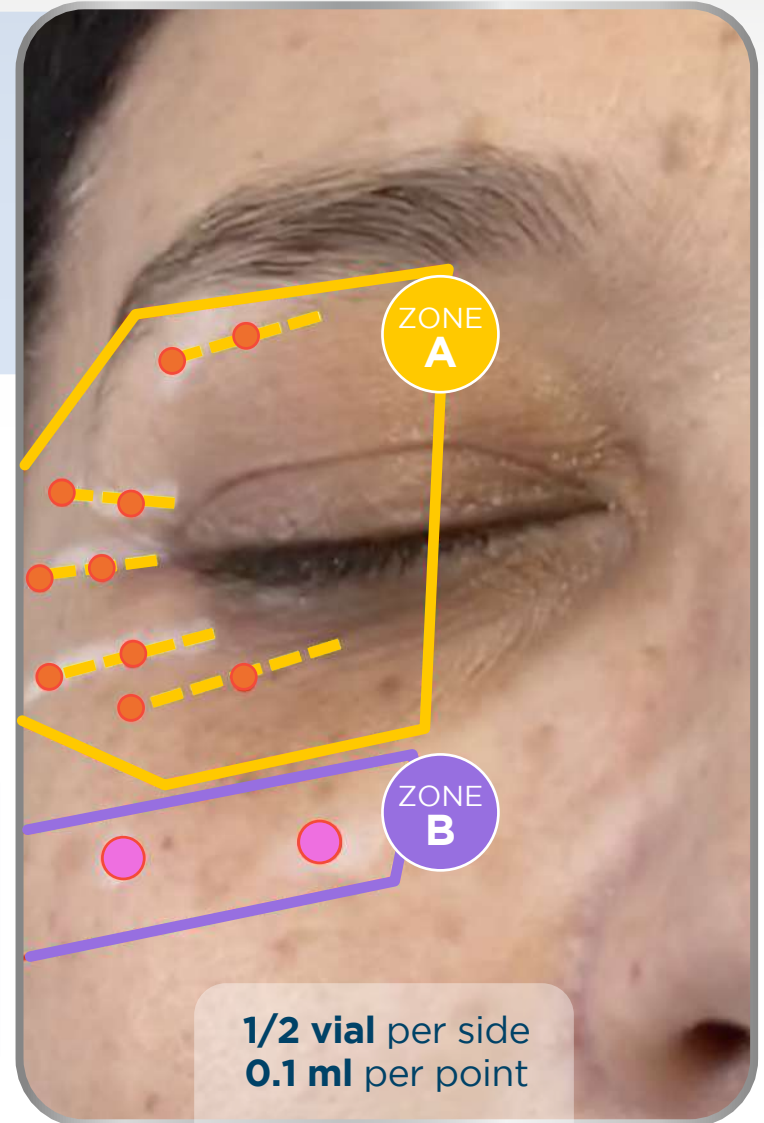

Sunekos[®] 200

Technique for Dark Eye Circles & Bags



BEFORE



**48H AFTER
TREATMENT**

Sunekos[®] 200

Dark Eye Circles



BEFORE



**1 WEEK AFTER
FIRST SESSION**

Sunekos® 200

Eye Bags



BEFORE



AFTER

Sunekos[®] 200

Eye Bags and Dark Eye Circles



BEFORE




AFTER

Sunekos[®] 200 + Sunekos[®] 1200 Hands

Thin skin of hands is one of the most affected by the loss of elastic properties.

The combined action of **Sunekos[®] 1200 + Sunekos[®] 200** injected with **Cushion Technique** helps to **rejuvenate the skin**.

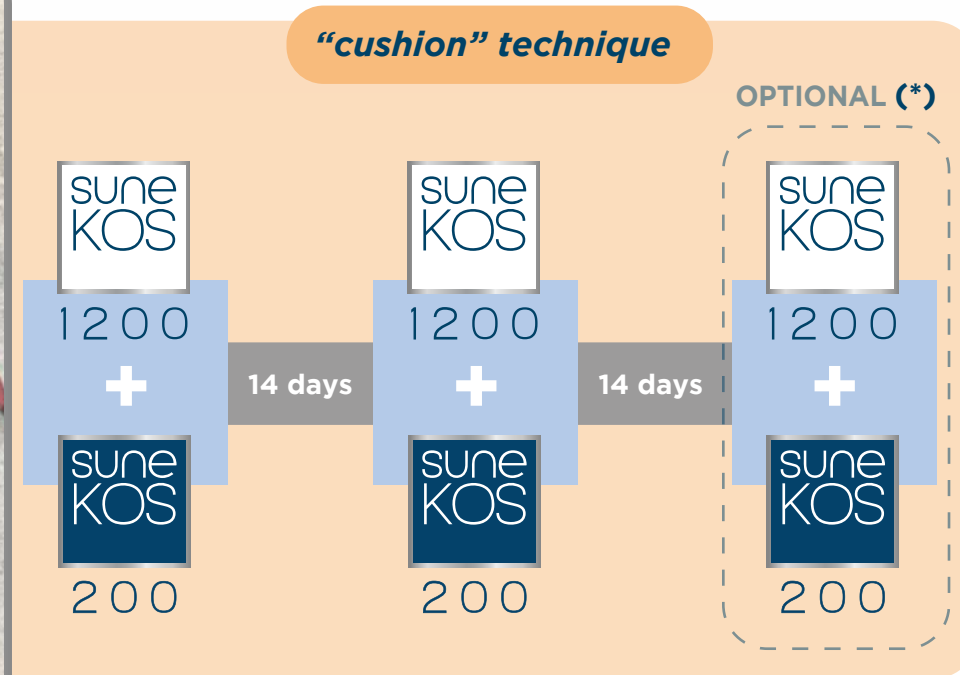

IT IS
RECOMMENDED
TO ADD
Sunekos Trophic Peeling
AS A BIO-RESTRUCTURING
SYSTEM BEFORE STARTING
THE PROTOCOL
AND AFTER THE LAST
INJECTION.



(*) to be evaluated by doctor but recommended

Sunekos[®] 200 + Sunekos[®] 1200

Cushion Technique for hands

- Cannula 25 or 27G and 50 mm
- Fan Technique

STEP
1

sune
KOS
1200



STEP
2

sune
KOS
200



- Needle 27-30G and 12-13 mm
- 0.1 ml injection at 1.5 cm distance
- Apply about 15 bolus avoiding veins

Sunekos® 200

Forehead

- Cannula 25 or 27G and 50 mm
- Deep Dermis injection
- Fan Technique
- ½ vial per side



The injection of **Sunekos® 200** with Cannula helps in reducing significantly forehead wrinkles in a natural way.

The procedure consists in 3 phases:

PHASE
1

Mechanical action of the Cannula to detach the tissue

PHASE
2

Injection of the product with the *FAN TECHNIQUE*

PHASE
3

Massage the treated area to stimulate the diffusion of the product

PHASE
1

IS PARTICULARLY IMPORTANT ESPECIALLY AT THE FIRST APPLICATION.

Sunekos® 200 counteract atrophy and tissue degeneration.

It can be used combined with other treatments improving their lasting and reducing the number of their applications.

sune
KOS
200

14 days

sune
KOS
200

14 days

sune
KOS
200

Sunekos[®] 200

Forehead



BEFORE



**1 WEEK AFTER
FIRST SESSION**

Sunekos[®] 200

Full Face



BEFORE



AFTER